

Applications Laboratory Thermophysical Properties Section

Thermal Conductivity of
Lewco Super Mat 2

6753-P-18_4



J. Hanss

August 22nd, 2018

NETZSCH-Gerätebau GmbH, Wittelsbacherstr. 42, 95100 Selb, Germany

Tel.: +49 9287 881 183, e-mail: ngb_laboratory@netzsch.com

Introduction

The Thermophysical Properties Section of the NETZSCH Applications Laboratory received Pyrogel XTE Brook from Lewco Speciality Products Inc., USA, for thermal conductivity determination.

The guarded hot plate system GPH 456 *Titan*[®] from NETZSCH is capable of operation between -160°C and 600°C (T_{mean} = sample temperature). The system is ideally suited for the measurement of insulating materials with thermal conductivities of up to 1 W/(m·K). The system is vacuum-tight by design, allowing measurements to be carried out under inert and oxidizing atmospheres as well as under vacuum. Thanks to the innovative plate material and plate design as well as the innovative temperature measurement, the system allows measurements with outstanding precision and repeatability. The plate stack is positioned in an environmental chamber avoiding radial heat flow. Furthermore, the environmental chamber compensates for heat flow via the wiring. Therefore, the system is a true absolute measurement tool for the characterization of insulating materials. The system allows measurements at well-defined temperature differences (ΔT) of 2 K to 40 K over the sample. The system runs fully automatically. No user input is required during measurement. The system works in accordance with ASTM C 177, ISO 8302, DIN EN 12667 and DIN EN 12939. System control is realized with state-of-the-art, 32-bit MS[®] Windows[™] software.



Further information about the GPH 456 *Titan*[®] can be found on the NETZSCH homepage:

<https://www.netzsch-thermal-analysis.com/en/products-solutions/thermal-diffusivity-conductivity/gph-456-titan/>

Experimental

The material tested was "Lewco Super Mat 2". The material was measured with spacers (thickness 15,3 mm) and minimum compression. Two specimens with about 300 mm x 300 mm were used for the test. The density during the measurement was approx. 181 kg/m³.

The measurements were performed under Nitrogen atmosphere (200 ml/min). The measured temperatures (T_{mean}) were 0°C, 100 °C, 200 °C, 300 °C and 600 °C (respectively 32 °F, 212 °F, 392 °F, 572 °F and 1112 °F) with a temperature difference across the specimens (ΔT) of 20 K except 1112°F, at this temperature only a ΔT of 10 K was possible.

Results

Figure 1 depicts the thermal conductivity of Lewco Super Mat 2 between 32 °F and 1112 °F. The thermal conductivity increases with increasing temperature as expected for a porous insulation material due to higher heat transfer by radiation at higher temperatures. Table 1 summarizes the corresponding measurement data.

After the measurement a mass loss of approx. 8 % was observed, which is most probably caused by decomposition of PTFE. Pyrolytic soot was found on the plates of the instruments and the color of the sample changed to dark grey/black after the measurement.

Table 1 summarizes the corresponding measurement data.

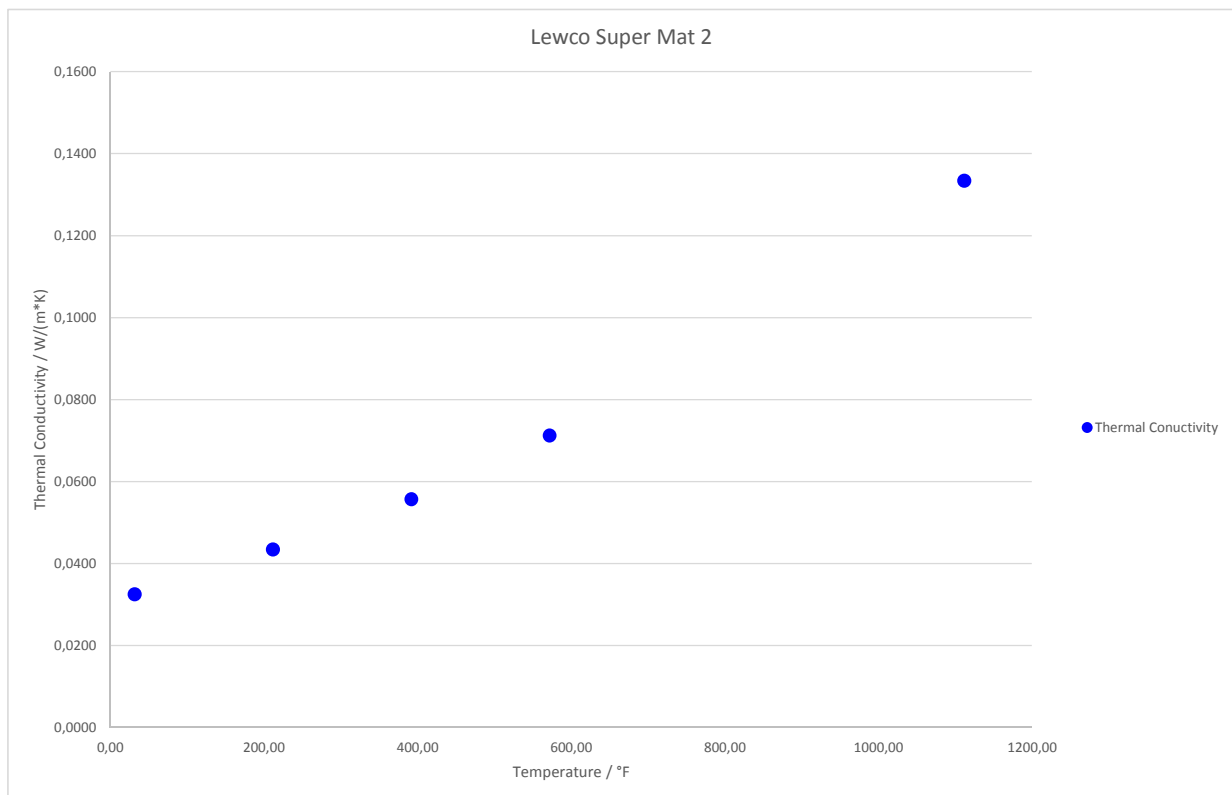


Figure 1: Thermal conductivity of Lewco Super Mat 2

Table 1: Thermal conductivity of Lewco Super Mat 2

Density @ RT ~ 181 kg/m³

Spacers (thickness) = 15.31 mm

$T_{\text{mean}} / ^\circ\text{C}$	$T_{\text{mean}} / ^\circ\text{F}$	$\Delta T / \text{K}$	Thermal Conductivity / $\text{W}/(\text{m}\cdot\text{K})$	Thermal Resistance / $(\text{m}^2\cdot\text{K})/\text{W}$
-0.02	31.97	20.04	0.0325	0.4706
100.02	212.04	19.96	0.0434	0.3524
200.02	392.04	19.95	0.0557	0.2747
300.02	572.04	19.95	0.0712	0.2149
600.02	1112.04	10.00	0.1334	0.1148

August 23rd, 2018
NETZSCH-Gerätebau GmbH
Applications Laboratory
i.V.



Jan Hanss